

PRESTO W91ttx

Powerful Temperature Control Systems for demanding applications with a temperature range from -91 °C to +250 °C

The PRESTO® W91 Highly Dynamic Temperature Control Systems are amongst the most powerful units available. Their impressive power is harnessed by cutting-edge control technology. The W91 provides extraordinary heating, cooling, and pumping performance. But it also comes with the functions and features that make the new generation of PRESTO® so unique.

Your advantages

- Extremely fast cool-down and heat-up times
- Wide working temperature range without changing the bath fluid
- Very fast compensation of exothermic and endothermic reactions
- Unsurpassed power and efficiency
- Heating capacity up to 36 kW
- Ambient temperature range +5 °C to +40 °C
- Space optimized design creates more room directly next to the units
- Integrated 5.7" industrial color touchscreen displays all essential information and enables simple fingertip control
- Extensive warning, protection, and monitoring functions with detailed self-explanatory messages
- ICC cascade control for extraordinary precision, temperature stability ± 0.05 °C ... ± 0.2 °C
- Integrated programmer with real-time clock
- Filling level indicator and pump capacity displayed electronically
- Powerful circulating pumps, electronically adjustable in stages or by setting the pressure value
- Interface for SD memory card
- Connections for USB, Ethernet, RS232, and Alarm output
- Optional analog connections, RS485, Profibus DP, Modbus



Made in Germany

Technical Data


| | | | | | | | | | |
|-----------------------------------|---|-----|-----|----|----|-----|-----|-----|-----|
| Order No. | 9421913.TT | | | | | | | | |
| Category | Temperature Control PRESTO | | | | | | | | |
| Working temperature range (°C) | -91 ... +250 | | | | | | | | |
| Temperature stability (°C) | 0.05 ... ± 0.2 | | | | | | | | |
| Setting / display resolution | 0.01 °C | | | | | | | | |
| Temperature Display | TFT Touchscreen | | | | | | | | |
| Heating capacity (kW) | 36 | | | | | | | | |
| Cooling capacity (Medium Ethanol) | °C | 200 | 100 | 20 | 0 | -20 | -40 | -60 | -80 |
| | kW | 11 | 11 | 11 | 10 | 9.5 | 9 | 6.5 | 1.5 |
| Pump capacity flow rate (l/min) | 18 ... 70 | | | | | | | | |
| Pump capacity flow pressure (psi) | 11.6...79.77 | | | | | | | | |
| Refrigerant | R23, R404A | | | | | | | | |
| External Pt100 sensor connection | integrated | | | | | | | | |
| Digital interface | RS232, RS485, Optional Profibus, SD memory card, USB, | | | | | | | | |

| | |
|--|--|
| | Ethernet, Modbus |
| Ambient temperature | 5 ... 40 °C |
| Dimensions W x L x H (inch) | 37.4 x 50 x 74.8 |
| Weight (LBS) | 1775 |
| Classification according to DIN12876-1 | Classification III (FL) according to DIN 12876-1 |
| Analog connection input / output | Optional |
| Cooling of compressor | 2-stage Water |
| Power requirement V / Hz / A | 3x 480/60/60 |
| Available voltage versions | 400V/3PNPE/50Hz 480V/3PPE/60Hz |


All data refers to the nominal voltage of 400 V, 3-phase, nominal frequency of 50 Hz and ambient temperature of +20 °C. Cooling capacity measured at max. pump stage. All pump data refers to a bath fluid with a specific density of 1 kg/dm³.


Characteristics

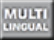
Display


 **State-of-the-art display technology**
TFT Display for comfortable user guidance, colored display of measurement values, graphs and control options, user-defined views

Operation


 **Optimal ease of use**
Touch screen for direct operation via display


 **Instructions inside**
Help menus and explanations in plain text for all control options, help messages and warning messages


 **Multilingual user guidance**
Language selection for display of control options, notifications and warning messages via touchscreen


 **Convenience for several users**
Administrator level for customizing instrument settings, user levels with limited permissions for fast and safe defined access, password protection, all levels adjustable

Temperature Control


 **For perfect results**
'Intelligent Cascade Control', automatic & self optimizing adjustment of PID control parameters, temperature stability ± 0.01 °C ... $\leq \pm 0.2$ °C


 **Full control**
'Temperature Control Features', for individual optimization, access to all important control parameters, additional settings for band limit, limits, co-speedfactor etc.


 **Control from the external application**
External Pt100 sensor connection for precise measurement and control directly in the external application

 **Highest measuring accuracy**
'Absolute Temperature Calibration' for manual compensation of a temperature difference, 3-point calibration


Refrigeration Technology


 **Consistent cooling capacity**
Easily removable venting grid for quick and easy cleaning


 **100 % Cooling capacity**
'Active Cooling Control' for cooling available throughout the entire working temperature range, fast cool-down even at higher temperatures


 **Energy saving cooling**
Proportional cooling control for automatic adjustment of cooling power or temporary switch-off of compressor as needed to save up to 90 % energy in comparison to unregulated cooling machines


Technical Features


 **Intelligent pump system**
Reliable and consistent pump capacity, electronically adjustable pump stages or pressure value, automatic adjustment of pump capacity to viscosity

 **Communication via networks**
For the remote control of instruments via Ethernet networks, full access to all functions of the unit via a networkcapable PC

 **Intelligent communication**
USB connection for data exchange (e.g. service data) or for wireless remote control via WirelessTEMP®

 **Data exchange via SD-Card**
For data exchange (e.g. service data) via SD memory card

 **Connections according to standard**
RS232/RS485 dual-interface for serial data transmission according to EIA-485 industry standard (2-wire bus technology), upgradable with Profibus DP

 **Comfortable program control**
Integrated programmer for the execution of time and temperature dependant profiles, 8 temperature profiles with 60 steps max., with real time clock



Quiet as a whisper

Efficient components produce only a minimal sound decibel level



Space-saving footprint

All connections as well supply and exhaust air are located at the front or rear, no venting grids on the sides, units can be placed close to each other or the application



Continuous operation up to +40 °C

Robust temperature control instrument, continuous operation even at ambient temperatures of up to +40 °C



Easy transport by one person

Ergonomic design facilitates moving and positioning by one person



Filling level at a glance

Backlit indicator for selected pump stages and filling volume



Quick support

If an error occurs, the integrated Black-Box function permits fast diagnosis by the JULABO service team

Warning & Safety Functions



Early warning system for high/low temperature limits

Maximum safety for applications, optical and audible signal when limits are exceeded.



Duplicate safety

Adjustable high temperature cut-off for internal tank and for integrated expansion vessel



For flammable bath fluid

Classification III (FL) according to DIN 12876-1

JULABO Contact

JULABO USA Inc.
884 Marcon Boulevard, Allentown, PA 18109
Phone: 610-231-0250
Phone: 800-458-5226 (Toll-Free)

info.us@julabo.com

JULABO Services

Product finder, Accessory search, 1 PLUS Warranty, Catalog download, Callback service, Operating manuals, Safety data sheets, Software and more

www.julabo.com